Specifications				
Specifications	Model No.			
220 V/230 V/240 V, 50 Hz (CE)	MDF-137-PE	MDF-237-PE	MDF-437-PE	
External dimensions (W x D x H)*1	640 x 687 x 881 (mm) 25.2 x 27.0 x 34.7 (inch)	905 x 687 x 881 (mm) 35.6 x 27.0 x 34.7 (inch)	1265 x 807 x 902 (mm) 49.8 x 31.8 x 35.5 (inch)	
Internal dimensions (W x D x H)	525 x 440 x 715 (mm) 20.7 x 17.3 x 28.1 (inch)	790 x 440 x 715 (mm) 31.1 x 17.3 x 28.1 (inch)	1140 x 550 x 733 (mm) 44.9 x 21.7 x 28.9 (inch)	
Capacity	138 liters (4.9 cu. ft.)	221 liters (7.8 cu. ft.)	425 liters (15.0 cu. ft.)	
Net Weight (Approx.)	52 <b>kg</b> (115 lbs.)	60 <b>kg</b> (132 lbs.)	81 <b>kg</b> (179 lbs.)	
Cooling performance	-30°C (AT: 35°C, no load)			
Control range	-20°C to -30°C			
Baskets	2	3	4	
Access ports	17mm diameter, right-hand side and bottom left			
Compressor	Hermetic, 150 W	Hermetic, 200 W	Hermetic, 250 W	
Refrigerant	HFC HFC and HFO combined			
Alarm system	• High/low temperature alarm (SV ±5°C to ±15°C adjustable) • Power failure • Remote alarm contact			
Calibration	Zero adjustment via control panel			
Options	Temperature recorder (MTR-G85C-PE, MTR-4015LH-PE)			
Baskets	MDF-13B2-PW/13B3-PW		MDF-43B2-PW/43B3-PW	

	Model No.		
110 V, 60 Hz	MDF-U5312-PT	MDF-U334-PT	_
220 V, 50 Hz	MDF-U5312-PB	_	MDF-U443-PB
220 V, 60 Hz	MDF-U5312-PK / MDF-U5312-PR	MDF-U334-PK/MDF-U334-PR	MDF-U443-PK/MDF-U443-PR
230 V/240 V, 50 Hz (CE)	MDF-U5312-PE	_	MDF-U443-PE
220 V/230 V/240 V, 50 Hz (CE)	_	MDF-U334-PE	_
External dimensions (W x D x H)*1	804 x 772 x 1802 (mm) 31.7 x 30.4 x 70.9 (inch)	614 x 709 x 1620 (mm) 24.2 x 27.9 x 63.8 (inch)	800 x 832 x 1810 (mm) 31.5 x 32.8 x 71.3 (inch)
Internal dimensions (W x D x H)	658 x 607 x 1272 (mm) 25.9 x 23.9 x 50.1 (inch)	490 x 486 x 1290 (mm) 19.3 x 19.1 x 50.8 (inch)	640 x 615 x 1090 (mm) 25.2 x 24.2 x 42.9 (inch)
Capacity	482 liters (17.0 cu. ft.)	274 liters (9.7 cu. ft.)	426 liters (15.0 cu. ft.)
Net Weight (Approx.)	134 <b>kg</b> (295 lbs.)	81 <b>kg</b> (179 lbs.)	213 <b>kg</b> (470 lbs.)
Cooling performance	-30°C (AT: 35°C, no load)		-40°C (AT: 35°C, no load)
Control range	-20°C to -30°C		-15°C to -40°C
Baskets	Top: 6 small, Bottom: 4 large	4 large, 1 small	_
Access ports	30mm diameter, rear	30mm diameter, left-hand side	40mm diameter, left-hand side
Compressor	Hermetic, 400 W (PT, PB)	Hermetic, 200 W	Hermetic, 450 W/750 W
Refrigerant	HFC		HFC and HFO combined
Alarm system	•High/low temperature alarm (SV ±5°C to ±15°C adjustable) •Power failure •Remote alarm contact •Part replacement notification	High/low temperature alarm (SV ±5°C to ±15°C adjustable)    Power failure    Part replacement notification    Remote alarm contact    Self diagnostics    Remote alarm contact	High/low temperature alarm (SV ±5°C to ±20°C adjustable) •Power failure •Filter clogging •Remote alarm contact •Battery/fan motor replacement notification
Calibration	Zero adjustment via control panel		_
Options	Temperature recorder (MTR-G85C-PE/MTR-G85A-PA, MTR-4015LH-PE/PA)		Temperature recorder (MTR-85H-PW, MTR-4015LH-PE/PA)

<sup>\*</sup>Cooling performance is indicated by the temperature reached at the center of the freezer (at ambient temperature of 30°C with no load). In order to use the freezer at a stable temperature for a long time, it is recommended that the temperature be set to at least 5°C higher than the indicated lowest temperature. In addition, depending on the usage conditions, it may not be possible to reach the indicated lowest temperature.

Caution: PHC Corporation guarantees this product under certain warranty conditions. However, please note that PHC Corporation shall not be responsible for any loss or damage to the contents of the

Appearance and specifications are subject to change without notice. \*1 External dimensions of main cabinet only - see dimension drawings showing handles and other external projections.

## **Options**

## Storage Baskets



For MDF-137/237 MDF-13B2-PW (set of two) MDF-13B3-PW (set of three) 255(W) x 420(D) x 205(H)mm






MDF-43B2-PW (set of two) MDF-43B3-PW (set of three) Outer dimensions 275(W) x 530(D) x 205(H)mm dimensions: 230(W) x 480(D) x 200(H)mm

## Temperature Recorders




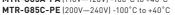



- Recorder fixing for MDF-137/237/437: MDF-S740-PW Recorder fixing for MDF-U334 upper room: MDF-S740T-PW
- Recorder fixing for MDF-U5312 lower room: MPR-S7-PW
- Recording paper: RP-G85-W Ink pen: PG-R-PW



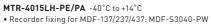

- Recorder fixing for MDF-U5312 lower room (no kit required for MDF-U334/MDF-U443]: MPR-S30-PW
- Recording paper: RP-40-PW





• Recording paper: RP-85-PW • Ink pen: DF-38FP-PW

Interface board (For data acquisition system MTR-5000 users only)

MTR-L03-PW / MTR-480-PW

# Preservation (freezers, refrigerators) and Culturing (incubators)



PHC Corporation, Biomedical Division

1-1-1 Sakada, Oizumi-machi, Ora-gun, Gunma 370-0596, Japan





PHC Corporation Biomedical Division is certified for: Environmental management system: IS014001

DISTRIBUTED BY:



https://www.phchd.com/global/biomedical/

Printed in Japan 1002-2018-04-CC

## **Biomedical Freezers**



# MDF-U443 無霜冷凍櫃



# **Biomedical Freezers**

## **Chest type**

## Power to cope with frequent door openings

## Small size suitable for personal use


- Digital temperature display as standard feature
- Compact, personal size



MDF-137 138 liters

## **Performance Data**

--- Pull-down characteristics under no-load condition Pull-up characteristics under no-load condition

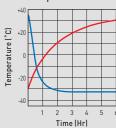

## **Features Enhanced Operation**

• Front-mounted display/control panel (with microprocessor) located at convenient waist height



- Memory backup
- Temperature display
- Front access calibration for temperature
- The control panel, alarm system and non-volatile memory are the same for all models in the series.

### Versatile mid-size model

- Digital temperature display as standard feature
- Easy-to-use 200-liter mid-size model



221 liters (7.8 cu.ft.

### Large size 400-liter class

- Digital temperature display as standard feature
- Large-capacity, economical size



425 liters

## • Single electrical box makes servicing simpler.

- The chamber temperature is displayed for five seconds if BUZZER key is depressed during power failure alarm.
- After a power outage, operation resumes at pre-outage settings (non-volatile memory for temperature and alarm temperature settings)
- Control panel can be reset to zero for
- Two access ports (right-hand side and bottom left) (chest type)
- Four casters and two adjustable feet

## Recording Features (Optional)

Recorder can be fitted in convenient position that will not interfere with installation. Choice of optional recorders

- 7-day circular recorder (MTR-G85C-PE) and 32-day strip recorder (MTR-4015LH-PE) are available. (except MDF-U443)
- Front access calibration for temperature recorder

### **Environment-Friendly**

- HFC refrigerants (non-CFC, non-HCFC, HFC and HFO combined) is environment-friendly but offers superb performance.
- Pre-coated metal body causes less environmental damage than conventional paint coating.

## **Upright type**

## Maximum capacity in minimum

### Large-capacity, economical, space-saving model

- Top/bottom twin door design prevents cold air escaping
- Convenient storage drawers



MDF-U5312

482 liters | 17.0 cu.ft.

## Easy-to-use, slim & compact design

- Microprocessor digital temperature display
- Convenient storage drawers



MDF-U334

274 liters

## Powerful fan ensures temperature uniformity throughout the freezer cabinet.

### Forced air circulation & auto defrost


- High freezing power increases temperature
- Auto defrost for easy maintenance



426 liters


## --- Pull-down characteristics

Pull-up characteristics during power failure

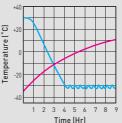

## Safety Features

- Alarm & Recording System
- High / low temperature alarm
- Power failure alarm
- Remote alarm contact
- Self diagnostics function. Should an abnormality be picked up, an error code will be displayed in red figures.
- Breaker switch turns power off in event of over current abnormality.
- Balance hinges keep the lid open in any position to make this chest freezer easier to use. This means you have both hands free to work. (chest type)

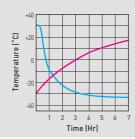

## • Secure lock safeguards the contents of the


freezer and keeps hazarous items out of harm's way.




Newly designed lid lock (chest type).



• Door latch can accommodate a pad lock for added protection of valuable

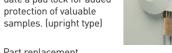


• Part replacement notification (MDF-U5312)

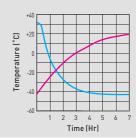

in cabinet Ambient temp. 35°C, No-load

## Hot line for secure sealing

Moisture condensation at the top edges of the cabinet due to differences in temperature inside and out causes frost and icing problems that may reduce heat insulation efficiency and obstruct door movements. These problems are prevented by the "hot line" by means of which hot gas from the higher temperature circuit is circulated through the problem areas.

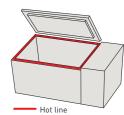

## **Dimensions**

